

## A. Operation Functions



### Display Keys

- 1. Power Key:** Turns the shoes on or off.
- 2. Massage Mode Key:** Allows free selection of 8 massage modes.
- 3. Intensity Key:** Adjusts the strength of the massage.
- 4. Language Key:** Selects the language mode (Chinese or English).
- 5. Time Duration Key:** Choose between 5, 10, or 15 minutes.

### Preparation:

- Clean your feet and remove any dirt from your skin.
- **Avoid wearing socks while using the shoes in massage mode**
- Use an alcohol wipe to **clean the shoe. Do not soak shoes in water**

### Usage:

- Connect the wire to the buttons on the shoe or massage patch.
- Start with the weakest pulse intensity level (*Level 1*) and gradually increase to *Level 15* to avoid sudden discomfort.
- Charge the product by inserting a USB into the charging hole.

## B. Safety Precautions

The company cannot be held responsible if you fail to use this product properly and safely in accordance with this provided instruction manual.



### 1. Safe Usage Guidelines:

- Do not use while driving.
- Remove any jewellery or items that may interfere with the massage area.
- Avoid using in humid environments such as bathrooms and spas.
- Keep the product away from water sources, fire, and corrosive environments.
- Always comply with safety regulations in your location.

### 2. Usage Restrictions:

- Avoid using the massage pads on areas of your body where you have undergone surgery.
- Do not apply to damaged, inflamed, or congested skin.
- **If you have heart disease, acute conditions, osteoporosis, or soft tissue diseases, refrain from using the product.**
- Avoid use if you experience discomfort or have slow movement, physical or mental discomfort, or mental retardation.
- Cease use immediately if you feel any discomfort.
- Not suitable for use during menstruation or pregnancy.
- Children under 6 years old should not use this instrument.



### 3. Cleaning Instructions:

- Clean the shoe using alcohol-soaked cotton; avoid water.
- Do not use sharp objects that could scratch or puncture the product.



## B. Precautions continued...

### 4. Product Usage:

- This product is a household electronic health device and is not intended for medical diagnosis, treatment, or commercial purposes.

### 5. Charging Safety:

- Charge the product only using the provided charging equipment.
- Be aware of electricity safety precautions during charging.

### 6. Abnormal Situations:

- If you notice any abnormalities, immediately disconnect the power supply.
- Do not attempt to disassemble or repair the product.

**Do not wear the shoes with the following medical electronic equipment**



- Heart rate adjusters and internal pulse rate medical electronic equipment
- Artificial cardiopulmonary pacemakers and life-sustaining medical electronic equipment
- Electrocardiogram equipment

**Note well:** Using this product alongside medical equipment (whether implanted or nearby) while operating the pulse massage product with the shoe simultaneously may endanger your body. It could affect the product's function and the medical equipment's performance.



## Caution:

Use with care if you: 

- Have a fever of 38°C or above.
- Have extremely sensitive skin.
- Are a heart disease patient (avoid using near the heart).
- Avoid placing the product near the mouth, eyes, perineum, or skin trauma areas.
- Prohibit using this product with other devices or applicators.

## Product Maintenance & Storage

Product Maintenance:

- Wipe the product with a soft cloth (avoid excessive wetness).
- Regularly check for wear, tear, and misuse.


Storage:

- Keep away from children, moisture, high temperatures, direct sunlight, and water splashes.



## Disposal & Recycling

- Switch off the external power supply.
- Unplug the power cord.
- Remove the battery.
- Dispose of the battery at local recycling points following environmental laws and guidelines.
- Remove the battery before disposing of the product.
- Do not attempt to repair or disassemble any parts.

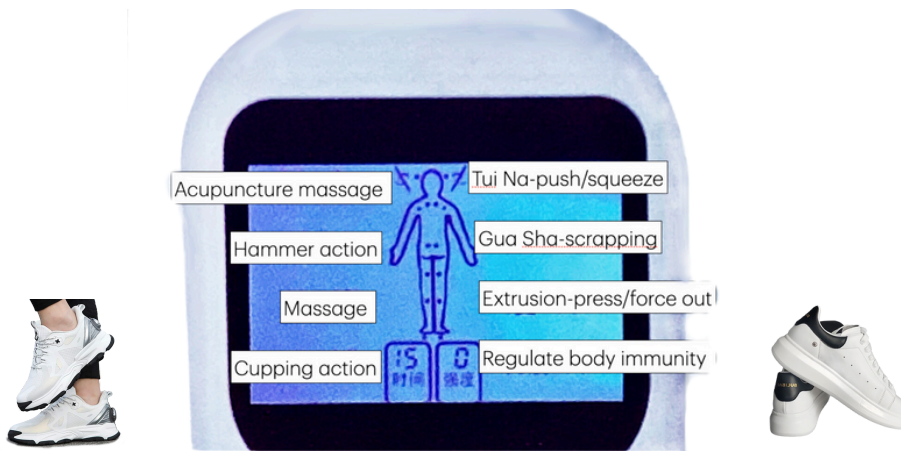

Product Name	Multifunctional Massage	Fault Status	Reason	Solution
Product Number	YHK-902-1	No display	Battery Exhausted	Recharge with USB Cable
Product Type	Rechargeable	Massage not activating	Button incorrectly inserted No direct skin contact	Insert and fasten properly
Voltage	3.7Vm	Weak pulse intensity of massage	Skin not in direct contact with sole or patch Dirty or oily skin or patch	Clean and dry your feet, <b>do not wear socks while massaging.</b> Clean your skin/patch
Power	4W	Massage paused	Battery exhausted	Recharge with USB cable
Battery	90mAh	Irritated, red itchy skin at massage areas	Allergic skin reaction	Do not use if you have skin allergies
Charging Time	3 Hrs			
Standards	GB4706.10-2008, GB4706.1-2005			
Power Input	100-240V AC			
Host Frequency	1-350Hz			
DC	5V 300mA			
Output Power	0.2W			
Product Size	12.8 x 6.2 x 1.64 cm			
Product Weight	0.19kg			